


Water Quality by Monitoring the Natural Organic Matter of Aquatic Systems



AN_P30; May 2016, Georgios Arnaoutakis, Anna Gakamsky, Dirk Näther

Introduction

Water in aquatic environments may consist of a complex mixture of organic compounds. This involves a continuum of natural organic matter of variable size, particulate or colloidal nature.^{1,2} The Natural Organic Matter (NOM) has been widely used to characterise water. Its measurement involves the Total Organic Content (TOC), the sum of particulate and Dissolved Organic Carbon (DOC),^{3,4} also known as humic substances or refractory organic substances.

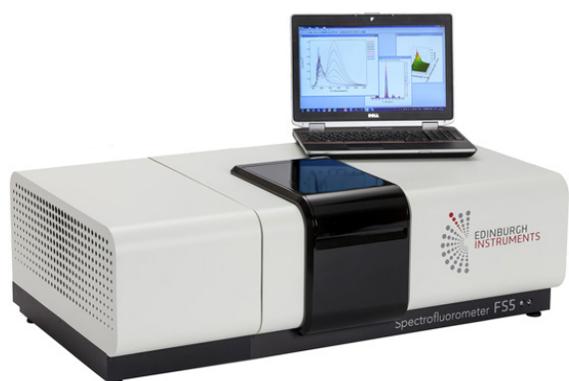

FSS Spectrofluorometer

Humic substances in aquatic systems originate from degradation of plant and animal tissue with a precursor being lignin which fluoresces at 360 nm upon excitation between 240 nm and 320 nm.⁵⁻⁷ Processes operative during the degradation are complex, with poly-condensation and formation of polyphenols being a possible candidate. The fluorescence can be related to substituted benzoic moieties with the start material being coniferyl alcohol, a basic unit of lignin.⁶ Model compounds and component analysis also associated this fluorescence to coniferyl alcohol, stilbene and phenyl-coumarone structures.⁷⁻¹⁰ On the other hand, protein fluorescence centres, especially observed in marine and pond water, are at the same wavelengths as those of tryptophan and tyrosine; although it is not known how these fluorescence centres relate to the structure of DOC.^{11,12} Fluorescence centres ascribed to humic-like and fulvic-like material occur at higher emission wavelengths.^{3,13}

The complexity of the absorption and emission spectra makes it impossible to clearly assign individual peaks occurring from independent chromophores. Different techniques are used to characterise NOM.⁴ A variety of analytical techniques such as fluorescence spectroscopy, Fourier Transform Infrared Spectroscopy, nuclear magnetic resonance, high pressure size exclusion chromatography, mass spectrometry have been widely employed,⁴ and even a combination of the techniques is used, especially when water components need to be quantified. In this application note, we present measurements of water obtained from river aquatic systems and show how fluorescence spectroscopy can easily provide initial results on the organic fingerprint of water.

Methods and Materials

Excitation and emission maps (EEM) were measured using a standard configuration FS5 Spectrofluorometer equipped with a 150 W standard xenon lamp and a photomultiplier tube detector (PMT-900). Higher diffraction orders were filtered by the integrated long wave-pass filters in the FS5. Excitation and emission bandwidth of 4 nm were used, while the integration time of 0.1 s and a step of 2 nm enabled the acquisition of a complete EEM in 60 min.

Acquisition times as low as 15 min were also obtained for the same parameters and an integration time of 1 ms.

Water samples were acquired from two locations upstream and downstream of the river Almond in Livingston, UK. A water sample from the tap was used as a reference. All three samples were filtered for particulates,¹⁴ in a syringe filter of 0.20 μm pore size (Sartorius, Minisart 17597) which resulted in OD=0.3 at 250 nm. The samples were measured in quartz cuvettes of 10 mm path-length in right angle geometry.

Results-Discussion

The EEM displayed in Figures 1, 2 and 3, correspond to a downstream, an upstream river sample and a water sample from the tap, respectively. The EEM have been scaled to the peak of visible DOC of the downstream river sample. All maps show the band of Raman scattering proportional to the excitation wavelength which can be further normalised and corrected for inner-filter to acquire pure DOC maps.^{6,15-17}

The assignment of EEM peaks is found in literature as the traditionally defined humic-like regions.^{4,11,19} These are labelled as A with $\lambda_{\text{exc}}/\lambda_{\text{em}}=260\text{ nm}/400\text{--}460\text{ nm}$, and as C with $\lambda_{\text{exc}}/\lambda_{\text{em}}=320\text{--}360\text{ nm}/420\text{--}460\text{ nm}$. Soil fulvic acid (D and E) with $\lambda_{\text{exc}}/\lambda_{\text{em}}=390\text{ nm}/509\text{ nm}$ and $\lambda_{\text{exc}}/\lambda_{\text{em}}=455\text{ nm}/521\text{ nm}$, phytoplankton (N) with $\lambda_{\text{exc}}/\lambda_{\text{em}}=280\text{ nm}/370\text{ nm}$ and tryptophan/protein-like (T) with $\lambda_{\text{exc}}/\lambda_{\text{em}}=275\text{ nm}/340\text{ nm}$ were not observed in these water samples.

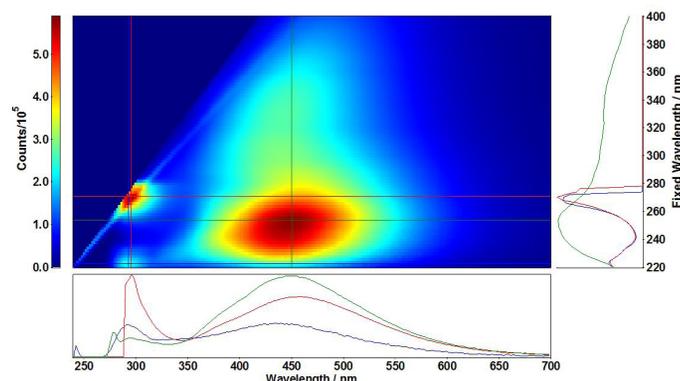


Figure 1: Excitation-map of water downstream of the river. The blue cross is at $\lambda_{\text{exc}}/\lambda_{\text{em}}=222\text{ nm}/292\text{ nm}$, the red at $\lambda_{\text{exc}}/\lambda_{\text{em}}=270\text{ nm}/296\text{ nm}$ and the green at $\lambda_{\text{exc}}/\lambda_{\text{em}}=254\text{ nm}/450\text{ nm}$

Water Quality by Monitoring the Natural Organic Matter of Aquatic Systems

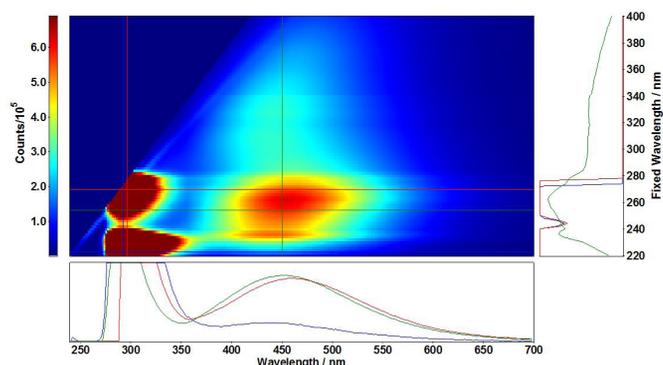


Figure 2: Excitation-emission map of water upstream of the river. Crossings are the same as shown in Figure 1.

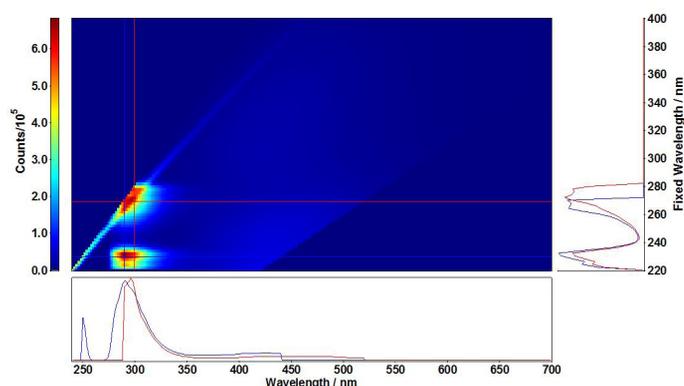


Figure 3: Excitation-emission map of tap water. The blue cross is at $\lambda_{exc}/\lambda_{em} = 230 \text{ nm}/290 \text{ nm}$, while the red at $\lambda_{exc}/\lambda_{em} = 270 \text{ nm}/300 \text{ nm}$

It can be seen that the DOC of water is centred at $\lambda_{exc}/\lambda_{em} = 250 \text{ nm}/450 \text{ nm}$, see green cross in Figures 1 and 2. This is in agreement with running river waters compared to the red-shifted EEM at $\lambda_{exc}/\lambda_{em} = 250 \text{ nm}/470 \text{ nm}$ from channels and wetlands.⁶

In addition to the humic-like peak, two additional peaks are present in all samples at $\lambda_{exc}/\lambda_{em} = 230 \text{ nm}/290 \text{ nm}$ and $\lambda_{exc}/\lambda_{em} = 270 \text{ nm}/300 \text{ nm}$, marked as blue and red in the figures, respectively. The peak with $\lambda_{exc}/\lambda_{em} = 230 \text{ nm}/290 \text{ nm}$ is associated with the superposition of electron-transfer and benzenoid bands in NOM molecules.²⁰

The UV humic-like peak with $\lambda_{exc}/\lambda_{em} = 270 \text{ nm}/300 \text{ nm}$ agrees very well with lignin units. In fact, lignin model compounds based on styrene derivatives with excitation 31847-36496 cm^{-1} (274-314 nm) and emission 28818-34246 \pm 200 cm^{-1} (292-347 nm) agree very well with the peaks in the presented maps.²¹

The broad peak at $\lambda_{exc}/\lambda_{em} = 320 \text{ nm}/440 \text{ nm}$ can also be associated to lignin. This agrees with the results of Radotic et al. on dehydrogenative polymers in water as another lignin model compound¹⁰, with excitation at 360 nm-465 nm and fluorescence at 360 nm-600 nm.¹⁰

The water supplied to domestic or industrial premises follows several treatment processes to make it usable or potable.¹³

Although the treatment process of the tap water is unknown to the authors, it can be seen that the humic-like peaks are absent in the EEM of tap water, indicating that certain DOM was removed. Similar fluorescence fingerprints have been reported for raw and treated water samples previously.¹³

Conclusion

It can be concluded from the presented excitation-emission maps that the content of aquatic systems is a complex mixture of degradation products of lignin containing aromatic, primarily benzene units. Moreover, it has been shown that the high sensitivity of an FS5 Spectrofluorometer enabled the resolution of both lignin, UV humic-like in addition to humic-like substances. The technique can also be used for rapid routine measurement of maps permitted by the high scanning speed of the instrument, as well as monitor the quality of water in treatment facilities.

References

- [1] Aiken, G. & Leenheer, J. Isolation and Chemical Characterization of Dissolved and Colloidal Organic Matter. *Chem. Ecol.* 8, 135–151 (1993).
- [2] Aiken, G. R., Gilmour, C. C., Krabbenhoft, D. P. & Orem, W. Dissolved Organic Matter in the Florida Everglades: Implications for Ecosystem Restoration. *Crit. Rev. Environ. Sci. Technol.* 41, 217–248 (2011).
- [3] Henderson, R. K. et al. Fluorescence as a potential monitoring tool for recycled water systems: A review. *Water Res.* 43, 863–881 (2009).
- [4] Matilainen, A. et al. An overview of the methods used in the characterisation of natural organic matter (NOM) in relation to drinking water treatment. *Chemosphere* 83, 1431–1442 (2011).
- [5] Kumke, M. U. & Frimmel, F. H. in *Refractory Organic Substances in the Environment* (eds. Frimmel, F. H. et al.) 215–231 (Wiley-VCH Verlag GmbH & Co. KGaA, 2002).
- [6] Hernes, P. J., Bergamaschi, B. A., Eckard, R. S. & Spencer, R. G. M. Fluorescence-based proxies for lignin in freshwater dissolved organic matter. *J. Geophys. Res. Biogeosciences* 114, G00F03 (2009).
- [7] Albinsson, B., Li, S., Lundquist, K. & Stomberg, R. The origin of lignin fluorescence. *J. Mol. Struct.* 508, 19–27 (1999).
- [8] S. INGEMAR FALKEHAG, JOSEPH MARTON & ERICH ADLER. in *Lignin Structure and Reactions* 59, 75–89 (AMERICAN CHEMICAL SOCIETY, 1966).
- [9] Lundquist, K., Josefsson, B. & Nyquist, G. Analysis of Lignin Products by Fluorescence Spectroscopy. *Holzforsch. - Int. J. Biol. Chem. Phys. Technol. Wood* 32, 27–32 (2009).
- [10] Radotic, K. et al. Component analysis of the fluorescence spectra of a lignin model compound. *J. Photochem. Photobiol. B* 83, 1–10 (2006).
- [11] Coble, P. G. Characterization of marine and terrestrial DOM in seawater using excitation-emission matrix spectroscopy. *Mar. Chem.* 51, 325–346 (1996).
- [12] Baker, A. et al. Measurement of protein-like fluorescence in river and waste water using a handheld spectrophotometer. *Water Res.* 38, 2934–2938 (2004).
- [13] Bagtho, S. A. Characterizing natural organic matter in drinking water treatment processes and trains. (UNESCO-IHE, Institute for Water Education, 2012).
- [14] Danielsson, L. G. On the use of filters for distinguishing between dissolved and particulate fractions in natural waters. *Water Res.* 16, 179–182 (1982).
- [15] Lawetz, A. J. & Stedmon, C. A. Fluorescence Intensity Calibration Using the Raman Scatter Peak of Water. *Appl. Spectrosc.* 63, 936–940 (2009).
- [16] Ohno, T. Fluorescence Inner-Filtering Correction for Determining the Humification Index of Dissolved Organic Matter. *Environ. Sci. Technol.* 36, 742–746 (2002)

Water Quality by Monitoring the Natural Organic Matter of Aquatic Systems



-
- [17] Spencer, R. G. M. et al. Discriminatory classification of natural and anthropogenic waters in two U.K. estuaries. *Sci. Total Environ.* 373, 305–323 (2007).
- [18] Singh, S. et al. Comparison of Two PARAFAC Models of Dissolved Organic Matter Fluorescence for a Mid-Atlantic Forested Watershed in the USA, Comparison of Two PARAFAC Models of Dissolved Organic Matter Fluorescence for a Mid-Atlantic Forested Watershed in the USA. *J. Ecosyst. J. Ecosyst.* 2013, 2013, e532424 (2013).the USA, Comparison of Two PARAFAC Models of Dissolved Organic Matter Fluorescence for a Mid-Atlantic Forested Watershed in the USA. . *Ecosyst. J. Ecosyst.* 2013, 2013, e532424 (2013).
- [19] Stedmon, C. A., Markager, S. & Bro, R. Tracing dissolved organic matter in aquatic environments using a new approach to fluorescence spectroscopy. *Mar. Chem.* 82, 239–254 (2003).
- [20] Korshin, G. V., Li, C.-W. & Benjamin, M. M. Monitoring the properties of natural organic matter through UV spectroscopy: A consistent theory. *Water Res.* 31, 1787–1795 (1997).
- [21] Korschin, H. et al. Fluorescence Characteristics of Lignin Model Compounds. I. Styrene Derivatives. *Acta Chem. Scand.* 30b, 262-266 (1976).

For more information, please contact:

T: +44 (0) 1506 425 300 E: sales@edinst.com

F: +44 (0) 1506 425 320 W: www.edinst.com